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Safety information.....	13
Regulatory information.....	16

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
1. Indicateurs de batterie gauche et droite
2. Smart Assistant*
3. Change program
4. Program in use
5. Program options
6. Edge Mode+: *Optimize sound quality on demand**
7. Programs
8. My Hearing: *Self Check, Find My Hearing Aids, additional sound settings*
9. Mute
10. Volume slider
11. Accessories: *View streaming options for available wireless accessories.*
12. Home
13. Health: *Monitor hearing and physical activity**
14. Learn: *Personalized tips, videos, and more*
15. More: *Account, About & privacy, TeleHear*

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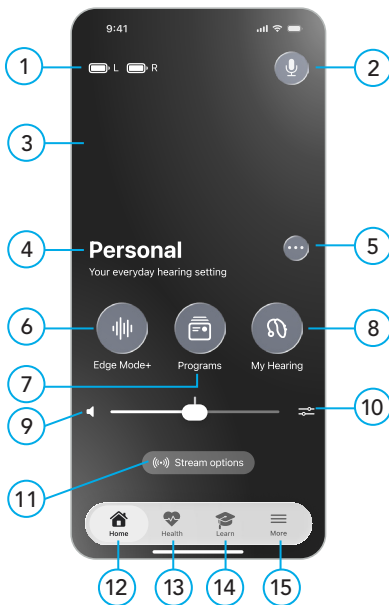1. On the **Home screen**, tap the **My Hearing button**.
2. The **My Hearing menu** appears. Tap Tinnitus.

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1. Ensure the **Volume tab** is open. (This tab should be open by default.)
2. Move the **right** (right hearing aid) and **left** (left hearing

aid) **sliders** to adjust the **volume** of the comfort sounds:

- For **LESS volume**, move sliders **LEFT**.
- For **MORE volume**, move sliders **RIGHT**.

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1. Appuyez sur le **onglet Vitesse**.
2. Move the **right** (right hearing aid) and **left** (left hearing aid) **sliders** to adjust the speed and variation of the comfort sounds:
 - For **LESS speed/variation**, move sliders **LEFT**.
 - For **MORE speed/variation**, move sliders **RIGHT**.

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Suivant :

After a brief countdown, this app sends an alert text message to your Fall Alert contacts. This contains details on your fall and a map of your location.

You'll also hear an **"alert sent" indicator** and receive a **confirmation message** that your alert was sent.

Your Fall Alert contact will receive an alert text message from this app, notifying them of your fall. After they receive it, you'll hear an **"alert received" indicator** through your hearing aids.

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1. On the **Home screen**, tap the **Health button**.
2. The **Health screen** appears. Tap **Wellness > Fall Alerts**.
3. Learn about Fall Alert on the first 2 **welcome screens**.
4. On the **third screen**, fill in your **contact's details**.
5. Tap **Send** to invite them to be your Fall Alert contact.

6. The **Invite Sent screen** confirms the invite message was texted to your contact. (It also allows you to add another Fall Alert contact.)

About & privacy

Your **About & privacy screen** lets you view info about this app, user guidelines and agreements. These include:

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It is good health practice for a person with a hearing loss to have a medical evaluation by a licensed physician (preferably a physician who specializes in diseases of the ear) before purchasing a prescription hearing aid. Licensed physicians who specialize in diseases of the ear are often referred to as otolaryngologists, otologists or otorhinolaryngologists. The purpose of the medical evaluation is to assure that all medically treatable conditions that may affect hearing are identified and treated before the hearing aid is purchased.

Following the medical evaluation, the physician will give you a written statement that states that your hearing loss has been medically evaluated and that you may be considered a candidate for a hearing aid. The physician will refer you to an audiologist or hearing aid dispenser, as appropriate, for a hearing aid evaluation.

The audiologist or hearing aid dispenser will conduct a hearing aid evaluation to assess your ability to hear with and without a hearing aid. The hearing aid evaluation will enable the audiologist or dispenser to select and fit a hearing aid to your individual needs.

If you have reservations about your ability to adapt to amplification, you should inquire about the availability of a trial-rental or purchase option program. Many hearing aid dispensers now offer programs that permit you to wear a hearing aid for a period of time for a nominal fee after which you may decide if you want to purchase the hearing aid.

In some geographies, you must have a medical evaluation before purchasing a prescription hearing aid. Some States allow an adult to waive the medical evaluation.

A hearing aid will not restore normal hearing and will not prevent or improve a hearing impairment resulting from organic conditions. Use of a hearing aid is only part of hearing habilitation and may need to be supplemented by auditory training and instruction in lip reading. In most cases infrequent use of a hearing aid does not permit a user to attain full benefit from it.

Some hearing instrument users have reported a buzzing sound in their hearing instrument when they are using mobile phones, indicating that the mobile phone and hearing instrument may not be compatible. It is well-known that mobile phones are potential sources of noise for hearing aids. Your Starkey Hearing Aids have been tested for compliance to two standards that define hearing aid immunity to digital wireless devices and meet the requirements of ANSI C63.19- 2019 as well as the criteria for user compatibility as defined by IEC 60118-13:2019.

 **CAUTION:** If set to the maximum output level and worn for periods of time exceeding the recommendations below, the patient's exposure to sound energy has the potential to exceed noise exposure limits. This instrument is intended for use for a maximum of sixteen (16) hours a day when set at the maximum output level.

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Multiflex Tinnitus Technology can be used as a part of a tinnitus treatment program.

Multiflex Tinnitus Technology plays a white noise through the hearing aid.

Multiflex Tinnitus Technology is programmed according to your hearing loss and preference, and your hearing professional can adjust the settings of Multiflex Tinnitus Technology to meet your needs.

Multiflex Tinnitus Technology may provide temporary relief of your tinnitus.

 **WARNING** Auto Alert may not detect 100 percent of falls

 **WARNING** Decreasing the Auto Alert Sensitivity may prevent some falls from being detected by your Fall Alert system. For example, Auto Alert may not detect a fall if:

- The Sensitivity setting is not appropriate for the user.
- The fall is very slow, or you slide down gradually.
- You get up and begin walking immediately after a fall.

As a reminder, you can initiate a Manual Alert if Auto Alert does not detect a fall. Manual Alert must be configured by your hearing professional before it can be used.

 **WARNING** Auto Alert may initiate false alerts. To prevent false-alert text messages from being sent to your contact(s), you may cancel the alert from either your smartphone or by pressing the user control on either hearing aid.

 **WARNING** To reduce Fall Alert communication failures:

- Your hearing aid(s) need to be powered on, paired and connected with your smartphone using Bluetooth
- The mobile device must be powered on, with the Hearing Aid app open (in the foreground or background).
- The mobile device must have a connection to the internet (via a cellular network or WiFi).

NOTE: Healthable features are intended for general wellness only. It is not intended to detect, diagnose, treat, cure, or prevent any disease or medical condition. Please consider the following:

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NOTE: For maximum comfort and reduced ear sensitivity, your hearing aids are intended to be worn up to 16 hours a day.

 **CAUTION:** The mobile application is software only and does not establish environmental operating conditions for the mobile device. The mobile device must be used in accordance with the manufacturer's specified environmental conditions.

 **CAUTION:** Use of the mobile application does not change hearing aid safety requirements. Follow all applicable hearing aid warnings and cautions.

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commerciales de Google LLC. Le mot-symbole Bluetooth® et les logos sont des marques déposées appartenant à Bluetooth SIG, Inc. et toute utilisation de ces marques par Starkey se fait sous licence. © m
©2026 Starkey Laboratories, Inc. Tous droits réservés. P00000145 8/26 BKLT3162-06-EN-XX