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Safety information.....	13
Regulatory information.....	16

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
1. Indicatori della batteria sinistra e destra
2. Smart Assistant*
3. Change program
4. Program in use
5. Program options
6. Edge Mode+: *Optimize sound quality on demand**
7. Programs
8. My Hearing: *Self Check, Find My Hearing Aids, additional sound settings*
9. Mute
10. Volume slider
11. Accessories: *View streaming options for available wireless accessories.*
12. Home
13. Health: *Monitor hearing and physical activity**
14. Learn: *Personalized tips, videos, and more*
15. More: *Account, About & privacy, TeleHear*

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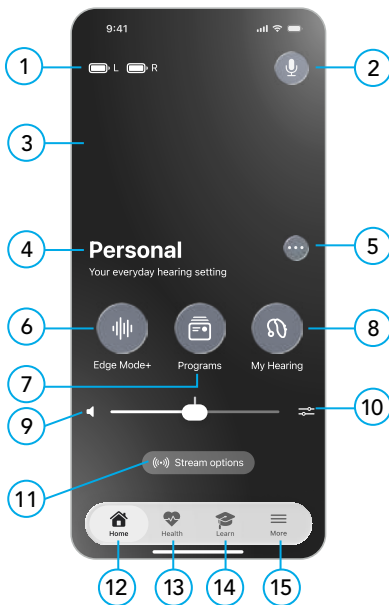
1. On the **Home screen**, tap the **My Hearing button**.
2. The **My Hearing menu** appears. Tap Tinnitus.

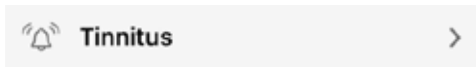3. Customize the **volume** or **speed** of the tinnitus relief sounds to suit your personal needs.

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1. Ensure the **Volume tab** is open. (This tab should be open by default.)
2. Move the **right** (right hearing aid) and **left** (left hearing aid) **sliders** to adjust the **volume** of the comfort sounds:

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Avanti:

After a brief countdown, this app sends an alert text message to your Fall Alert contacts. This contains details on your fall and a map of your location.

You'll also hear an **"alert sent" indicator** and receive a **confirmation message** that your alert was sent.

Your Fall Alert contact will receive an alert text message from this app, notifying them of your fall. After they receive it, you'll hear an **"alert received" indicator** through your hearing aids.

Canceling a Fall Alert

Cancel an **alert text message** to your contact by pushing **any button** on your hearing aids. Or tap **Cancel, I'm okay** on the **timer screen** before the countdown ends. When

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canceled, you'll hear an **"alert canceled" indicator** through your hearing aids.

Invite a contact

Set up Fall Alerts by inviting a contact who can be notified by an alert text message if you fall.

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1. On the **Home screen**, tap the **Health button**.
2. The **Health screen** appears. Tap **Wellness > Fall Alerts**.
3. Learn about Fall Alert on the first 2 **welcome screens**.
4. On the **third screen**, fill in your **contact's details**.
5. Tap **Send** to invite them to be your Fall Alert contact.
6. The **Invite Sent screen** confirms the invite message

was texted to your contact. (It also allows you to add another Fall Alert contact.)

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It is good health practice for a person with a hearing loss to have a medical evaluation by a licensed physician (preferably a physician who specializes in diseases of the ear) before purchasing a prescription hearing aid. Licensed physicians who specialize in diseases of the ear are often referred to as otolaryngologists, otologists or otorhinolaryngologists. The purpose of the medical evaluation is to assure that all medically treatable conditions that may affect hearing are identified and treated before the hearing aid is purchased.

Following the medical evaluation, the physician will give you a written statement that states that your hearing loss has been medically evaluated and that you may be considered a candidate for a hearing aid. The physician will refer you to an audiologist or hearing aid dispenser, as appropriate, for a hearing aid evaluation.

The audiologist or hearing aid dispenser will conduct a hearing aid evaluation to assess your ability to hear with and without a hearing aid. The hearing aid evaluation will enable the audiologist or dispenser to select and fit a hearing aid to your individual needs.

If you have reservations about your ability to adapt to amplification, you should inquire about the availability of a trial-rental or purchase option program. Many hearing aid dispensers now offer programs that permit you to wear a hearing aid for a period of time for a nominal fee after which you may decide if you want to purchase the hearing aid.

In some geographies, you must have a medical evaluation before purchasing a prescription hearing aid. Some States allow an adult to waive the medical evaluation.

A hearing aid will not restore normal hearing and will not prevent or improve a hearing impairment resulting from organic conditions. Use of a hearing aid is only part of hearing habilitation and may need to be supplemented by auditory training and instruction in lip reading. In most cases infrequent use of a hearing aid does not permit a user to attain full benefit from it.

Some hearing instrument users have reported a buzzing sound in their hearing instrument when they are using mobile phones, indicating that the mobile phone and hearing instrument may not be compatible. It is well-known that mobile phones are potential sources of noise for hearing aids. Your Starkey Hearing Aids have been tested for compliance to two standards that define hearing aid immunity to digital wireless devices and meet the requirements of ANSI C63.19- 2019 as well as the criteria for user compatibility as defined by IEC 60118-13:2019.

 **CAUTION:** If set to the maximum output level and worn for periods of time exceeding the recommendations below, the patient's exposure to sound energy has the potential to exceed noise exposure limits. This instrument is intended for use for a maximum of sixteen (16) hours a day when set at the maximum output level.

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As a reminder, you can initiate a Manual Alert if Auto Alert does not detect a fall. Manual Alert must be configured by your hearing professional before it can be used.

 **WARNING** Auto Alert may initiate false alerts. To prevent false-alert text messages from being sent to your contact(s), you may cancel the alert from either your smartphone or by pressing the user control on either hearing aid.

 **WARNING** To reduce Fall Alert communication failures:

- Your hearing aid(s) need to be powered on, paired and connected with your smartphone using Bluetooth
- The mobile device must be powered on, with the Hearing Aid app open (in the foreground or background).
- The mobile device must have a connection to the internet (via a cellular network or Wi-Fi).

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©2026 Starkey Laboratories, Inc. Tutti i diritti riservati. 8/26 BKL23162-06-IT-XX